



VetScalpel® CO₂ Laser
Pain management begins with less painful surgical procedures.

Since 1991, VetScalpel has provided veterinarians an affordable, Hollow Waveguide (HWG) CO₂ laser. We are the #1 choice of first-time buyers, and remain the #1 choice of second-time buyers.


VetScalpel is the 4th generation of our HWG CO₂ lasers. Our only focus is CO₂ (10,600 nm) lasers for soft tissue surgery.

- Enhanced surgical options
- Less painful procedures
- Ablation/vaporization tools never before available
- Practice differentiation



WHY CHOOSE VetScalpel?

- Less pain, better healing than scalpel (*studies upon request*)
- Limited hemorrhage improves visualization
- SuperPulse provides precision cutting
- Ergonomic HWG feels like a scalpel
- 4th generation, 35+ years experience

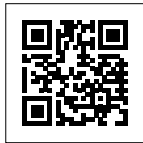

WATCH 20+ VIDEOS
www.vetscalpel.com/video




75+ CASE HISTORIES
www.vetscalpel.com/case-studies



Handpieces for every procedure

STEVE NIELSEN
 Int'l Business Development
sdnmt@gmail.com
 +1 404.626.3017

EVOLUTION OF CO₂ LASER SURGERY

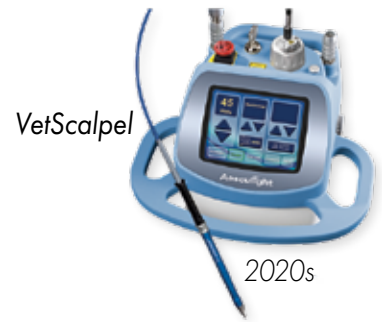

Why CO₂ Laser For Soft Tissue Surgery?

- Soft tissue is 60% - 85% H₂O
- Need wavelength with high H₂O absorption
- 10,600 nm has 1000 times better absorption than 810 or 980 nm in soft tissue
- Need wavelength with optimal, adjustable coagulation



VetScalpel® Features

Flexible Hollow Fiber Waveguide ensures surgical dexterity

20, 30, 40 Watt models with SuperPulse

Touch screen control panel ensures easy navigation

Continuous, Repeat Pulsing & Single Pulse modes provide precision and control

Unique SuperPulse mode enables tissue cooling while minimizing thermal damage

Surgeon-defined custom settings for fast recall

Autoclavable adjustable handpieces (0.25, 0.4 & 0.8 mm)

Fiber calibration port ensures power accuracy of the laser beam

Sealed metal laser tubes enable long life & low maintenance

Active forced-air cooling eliminates overheating

VetScalpel and plume evacuator are integrated for enhanced safety



UPCOMING SEMINARS/EVENTS

www.vetscalpel.com/events/



WHAT VETERINARIANS ARE SAYING

www.vetscalpel.com/testimonials/



REQUEST A QUOTE

www.vetscalpel.com/quote/

FOR VETERINARY USE ONLY

VETSCALPEL.COM
info@vetscalpel.com

Copyright ©2023 MTY. All Rights Reserved.



Bothwell, WA, USA